

*Contains Nonbinding Recommendations*

*Draft – Not for Implementation*

# Remanufacturing of Medical Devices

## Draft Guidance for Industry and Food and Drug Administration Staff

### *DRAFT GUIDANCE*

**This draft guidance document is being distributed for comment purposes only.**

**Document issued on June 24, 2021.**

You should submit comments and suggestions regarding this draft document within 90 days of publication in the *Federal Register* of the notice announcing the availability of the draft guidance. Submit electronic comments to <https://www.regulations.gov>. Submit written comments to the Dockets Management Staff, Food and Drug Administration, 5630 Fishers Lane, Room 1061, (HFA-305), Rockville, MD 20852. Identify all comments with the docket number listed in the notice of availability that publishes in the *Federal Register*.

For questions about this document regarding CDRH-regulated devices, contact the Regulation, Policy, and Guidance Staff at [RPG@fda.hhs.gov](mailto:RPG@fda.hhs.gov). For questions about this document regarding CBER-regulated devices, contact the Office of Communication, Outreach, and Development (OCOD) at 1-800-835-4709 or 240-402-8010, or by email at [ocod@fda.hhs.gov](mailto:ocod@fda.hhs.gov).



U.S. Department of Health and Human Services  
Food and Drug Administration  
Center for Devices and Radiological Health  
Center for Biologics Evaluation and Research

***Contains Nonbinding Recommendations***

***Draft – Not for Implementation***

# **Preface**

## **Additional Copies**

### **CDRH**

Additional copies are available from the Internet. You may also send an email request to [CDRH-Guidance@fda.hhs.gov](mailto:CDRH-Guidance@fda.hhs.gov) to receive a copy of the guidance. Please use the document number 17048 and complete title of the guidance in your request.

### **CBER**

Additional copies are available from the Center for Biologics Evaluation and Research (CBER), Office of Communication, Outreach, and Development (OCOD), 10903 New Hampshire Ave., Bldg. 71, Room 3128, Silver Spring, MD 20993-0002, or by calling 1-800-835-4709 or 240-402-8010, by email, [ocod@fda.hhs.gov](mailto:ocod@fda.hhs.gov) or from the Internet at <https://www.fda.gov/vaccines-blood-biologics/guidance-compliance-regulatory-information-biologics/biologics-guidances>.

***Contains Nonbinding Recommendations******Draft – Not for Implementation*****Table of Contents**

I. Introduction.....	1
II. Background.....	2
A. FDA activities.....	2
B. FDA’s current thinking .....	3
III. Scope .....	3
IV. Definitions .....	4
V. Guiding Principles.....	5
VI. Relevant considerations to determine if activities are remanufacturing .....	7
A. What is a significant change to device performance or safety specifications? .....	7
B. Determining whether activities are “remanufacturing” .....	8
A1. Add, remove, or change a component/part/material that directly or indirectly contacts body tissue?.....	12
A1.1 Is there a significant change to device performance or safety specifications?.....	12
A2. Add or remove component/part/material or change the dimensional or performance specifications of a component/part/material? .....	13
A2.1 Is there a significant change to device performance or safety specifications?.....	14
A3. Is there a new or modified risk or is there a change in the performance or safety specifications? .....	14
A3.1 Is there a significant change to device performance or safety specifications?.....	15
VII. Changes involving software.....	15
VIII.Considerations for labeling.....	16
Appendix A. Examples.....	18
(1) Component/part/material activities .....	18
(2) Software activities .....	28
Appendix B. Documentation examples.....	30

*Contains Nonbinding Recommendations*

*Draft – Not for Implementation*

# Remanufacturing of Medical Devices

## Draft Guidance for Industry and Food and Drug Administration Staff

*This draft guidance, when finalized, will represent the current thinking of the Food and Drug Administration (FDA or Agency) on this topic. It does not establish any rights for any person and is not binding on FDA or the public. You can use an alternative approach if it satisfies the requirements of the applicable statutes and regulations. To discuss an alternative approach, contact the FDA staff or Office responsible for this guidance as listed on the title page.*

### I. Introduction

Medical devices encompass a vast array of products with different technologies, product lifecycles, complexity, intended users, and environments of use. Many devices are reusable and need preventive maintenance and repair during their useful life. For these devices, proper servicing is critical to their continued safe and effective use. However, there is a lack of clarity regarding the distinction between “servicing” and “remanufacturing” of a device. Most notably, remanufacturing has implications for the regulatory responsibilities of entities performing these activities.<sup>1</sup>

This draft guidance is intended to help clarify whether activities performed on devices are likely “remanufacturing.” Such clarification is intended to help provide consistency and better understanding of applicable statutory and regulatory requirements. This draft guidance also includes recommendations for information that should be included in labeling to help assure the continued quality, safety, and effectiveness of devices that are intended to be serviced over their useful life. In drafting this guidance, FDA considered objective evidence and information learned from the Agency’s activities discussed in this draft guidance.

For the current edition of the FDA-recognized standard(s) referenced in this document, see the [FDA Recognized Consensus Standards Database](#).<sup>2</sup> For more information regarding use of consensus standards in regulatory submissions, please refer to the FDA guidance titled “[Appropriate Use of Voluntary Consensus Standards in Premarket Submissions for Medical](#)

<sup>1</sup> FDA’s [Report on the Quality, Safety, and Effectiveness of Servicing of Medical Devices](#) (FDA Report on Device Servicing) discusses medical device servicing in more detail. Available at <https://www.fda.gov/media/113431/download>.

<sup>2</sup> Available at <https://www.accessdata.fda.gov/scripts/cdrh/cfdocs/cfStandards/search.cfm>.

***Contains Nonbinding Recommendations******Draft – Not for Implementation***

Devices”<sup>3</sup> and “Standards Development and the Use of Standards in Regulatory Submissions Reviewed in the Center for Biologics Evaluation and Research.”<sup>4</sup>

The contents of this document do not have the force and effect of law and are not meant to bind the public in any way, unless specifically incorporated into a contract. This document is intended only to provide clarity to the public regarding existing requirements under the law. FDA guidance documents, including this guidance, should be viewed only as recommendations, unless specific regulatory or statutory requirements are cited. The use of the word should in Agency guidance means that something is suggested or recommended, but not required.

## **II. Background**

### **A. FDA activities**

FDA has been working to gain additional perspectives on the distinction between “servicing” and “remanufacturing” and has undertaken several efforts to help promote clarity. FDA opened a docket for public comment<sup>5</sup> and held a public workshop in 2016.<sup>6</sup> FDA received comments, complaints, and adverse event reports alleging inadequate servicing, which were discussed and analyzed in the FDA Report on Device Servicing,<sup>7</sup> published by FDA in May 2018 in accordance with section 710 of the FDA Reauthorization Act of 2017 (FDARA) (Pub. L. 115-52).<sup>8</sup>

In the FDA Report on Device Servicing, FDA concluded that a majority of the comments, complaints, and adverse event reports received by the Agency that referred to inadequate “servicing” causing or contributing to adverse events and deaths actually pertained to “remanufacturing.” This conclusion was based on FDA’s evaluation of the available objective evidence<sup>9</sup> related to the quality, safety, and effectiveness of medical device servicing.

In 2018, FDA released a white paper, opened a public docket, and held a public workshop to facilitate public discussion on the distinction between servicing and remanufacturing.<sup>10</sup> The

<sup>3</sup> Available at <https://www.fda.gov/regulatory-information/search-fda-guidance-documents/appropriate-use-voluntary-consensus-standards-premarket-submissions-medical-devices>.

<sup>4</sup> Available at <https://www.fda.gov/regulatory-information/search-fda-guidance-documents/standards-development-and-use-standards-regulatory-submissions-reviewed-center-biologics-evaluation>.

<sup>5</sup> 81 FR 11477. Public comments submitted to the docket are searchable under FDA-2016-N-0436, available at <https://www.regulations.gov/docket?D=FDA-2016-N-0436>.

<sup>6</sup> 81 FR 46694.

<sup>7</sup> Available at <https://www.fda.gov/media/113431/download>.

<sup>8</sup> FDA’s conclusions in this report were based on the available information, which included but was not limited to the information presented at the 2016 public workshop, responses to the docket request for comments, and evaluation of objective evidence related to the quality, safety, and effectiveness of medical device servicing.

<sup>9</sup> The objective evidence evaluated in the FDA Report on Device Servicing included a numerical estimation of service and repair entities, literature review, ECRI Institute analysis, medical device reports (MDR), and complaints that FDA received.

<sup>10</sup> Available at <https://wayback.archive-it.org/7993/20201222125933/https://www.fda.gov/medical-devices/workshops-conferences-medical-devices/public-workshop-medical-device-servicing-and-remanufacturing-activities-december-10-11-2018-12102018>. FDA requested comments through docket number FDA-2018-N-3741.

***Contains Nonbinding Recommendations***

***Draft – Not for Implementation***

white paper described FDA’s initial thoughts about guiding principles, provided a flowchart with accompanying text for understanding the distinctions, and contained a complementary approach for software, as well as considerations for labeling, and examples utilizing the flowchart. FDA also included targeted questions throughout the white paper for which the Agency sought feedback. FDA considered the comments from the public docket and discussions during the public workshop in developing this draft guidance.

## **B. FDA’s current thinking**

The distinction between “remanufacturing” and “servicing” is important to understand. Remanufacturing is the processing, conditioning, renovating, repackaging, restoring, or any other act done to a finished device that significantly changes the finished device’s performance or safety specifications, or intended use.<sup>11</sup> For the purposes of this guidance, we refer to the original equipment manufacturer’s (OEM’s) legally marketed finished device as the “legally marketed device.”

Servicing is the repair and/or preventive or routine maintenance of one or more parts in a finished device, after distribution, for purposes of returning it to the safety and performance specifications established by the OEM and to meet its original intended use.<sup>12</sup> As described in the FDA Report on Device Servicing, FDA’s authority to regulate the servicing of medical devices by any entity is grounded in the Agency’s authority to regulate medical devices and radiation-emitting electronic products under the Federal Food, Drug, and Cosmetic (FD&C) Act.

Irrespective of an entity’s self-identified designation as a “servicer” or “remanufacturer,” FDA focuses on the specific activities an entity performs on a particular device.<sup>13</sup> The determination of whether the activities an entity performs are remanufacturing affects the applicability and enforcement of regulatory requirements under the FD&C Act and its implementing regulations. FDA has consistently enforced requirements under the FD&C Act and its implementing regulations on entities engaged in remanufacturing, including but not limited to registration and listing, adverse event reporting, the Quality System (QS) regulation, and marketing submissions.

## **III. Scope**

Because of the apparent confusion between servicing and remanufacturing among entities performing these activities, FDA committed in the FDA Report on Device Servicing to issue guidance that clarifies the difference between servicing and remanufacturing activities. To assist

---

<sup>11</sup> See 21 CFR 820.3(w).

<sup>12</sup> For purposes of the report that Congress required FDA to post on its website, section 710(c) of FDARA (Pub. L. 115-52, 131 Stat. 1068) defines servicing to include, “with respect to a device, refurbishing, reconditioning, rebuilding, remarketing, repairing, remanufacturing, or other servicing of the device.” However, for purposes other than this report, FDA does not consider remanufacturing to be a type of servicing.

<sup>13</sup> The designations of servicer and remanufacturer are not mutually exclusive. The same entity may meet the definition of either designation based on their activities on one or multiple devices.

***Contains Nonbinding Recommendations******Draft – Not for Implementation***

with this clarification, FDA focuses this draft guidance on those activities that are likely remanufacturing.

This draft guidance addresses activities performed on devices that are intended to be reused and maintained. This draft guidance discusses whether activities performed by OEMs and third parties on such devices are likely remanufacturing. This draft guidance is not intended to adopt significant policy changes, but to clarify FDA’s current thinking on applicable definitions, and clarify, not change, the regulatory requirements applicable to remanufacturers. The concepts in this draft guidance are also not intended to alter or supersede existing regulations and policies related to the regulatory threshold for submitting a marketing submission for a device.

The products included within the scope of this guidance are devices as defined in section 201(h) of the FD&C Act, including software and electronic products that meet the definition of a device. In general, the concepts discussed in this guidance are meant to apply to all reusable devices, irrespective of their classification into class I, II, or III, including those subject to premarket approval. This guidance is not intended to address reprocessed single-use devices.

## **IV. Definitions**

The following definitions apply for the purposes of this guidance.<sup>14</sup>

- **Manufacturers (Manufacturers, OEMs, or Remanufacturers):** A manufacturer is any person who designs, manufactures, fabricates, assembles, or processes a finished device.<sup>15</sup> A remanufacturer is any person who processes, conditions, renovates, repackages, restores, or does any other act to a finished device that significantly changes the finished device’s performance or safety specifications, or intended use.<sup>16</sup> Remanufacturers are considered to be manufacturers.<sup>17</sup> For electronic products, a manufacturer is any person engaged in the business of manufacturing, assembling, or importing electronic products.<sup>18</sup>
- **Intended use:** The general purpose of the device or its function, which encompasses the indications for use.<sup>19</sup>
- **Performance specifications:** The performance characteristics of a device established by the OEM for the device to perform as intended, including those listed in device labeling or in finished product release specifications. Some examples include measurement accuracy, output accuracy, energy output level, and stability criteria.
- **Recondition/Refurbish/Rebuild:** Restores a medical device to the OEM’s original specifications or to be “like new.” The device may be brought to current specifications if the change(s) made to the device do not significantly change the finished device’s performance or safety specifications, or intended use. These

<sup>14</sup> Consistent with FDA’s current thinking in this context, some of the definitions that appeared in the FDA Report on Device Servicing have been modified to reflect updated understanding and practice.

<sup>15</sup> 21 CFR 820.3(o).

<sup>16</sup> 21 CFR 820.3(w).

<sup>17</sup> 21 CFR 820.3(o) and 820.3(w).

<sup>18</sup> 21 CFR 1000.3(n).

<sup>19</sup> FDA uses this term consistent with the meaning of intended uses in 21 CFR 801.4.



## *Contains Nonbinding Recommendations*

### *Draft – Not for Implementation*

- activities include repair of components, installation of OEM provided updates and upgrades, and replacement of worn parts.
- Remanufacture: Process, condition, renovate, repackage, restore, or any other act done to a finished device that significantly changes the finished device’s performance or safety specifications, or intended use.<sup>20</sup>
  - Repair: A type of servicing that returns a component to original specifications, including replacing non-working components or parts outside of routine or periodic upkeep for the current owner of the device.
  - Reprocessing: With respect to reusable devices, means validated processes used to render a medical device, which has been previously used or contaminated, fit for a subsequent single use on a patient. These processes are designed to remove soil and contaminants by cleaning and to inactivate microorganisms by disinfection or sterilization.<sup>21</sup>
  - Safety specifications: The safety characteristics of a device established by the OEM for the safe use of the device, including those incorporated into the device design and finished product release specifications, generally including the device’s compensating controls and risk mitigations. Some examples include alarms, sensors, and locking or fail-safe mechanisms.
  - Service: Repair and/or preventive or routine maintenance of one or more parts in a finished device, after distribution, for purposes of returning it to the safety and performance specifications established by the OEM and to meet its original intended use. Servicing excludes activities that significantly change the finished device’s safety or performance specifications, or intended use.
  - Third party servicers and Independent Service Organizations (ISOs): Entities, other than the OEM or healthcare delivery organizations, that maintain, restore, refurbish, or repair a finished device after distribution, for purposes of returning it to the safety and performance specifications established by the OEM and to meet its original intended use.

## **V. Guiding Principles**

In using this guidance to help determine whether activities are remanufacturing, FDA recommends application of the following guiding principles:

1. **Assess whether there is a change to the intended use** – Given that the purpose of servicing is to return the device to the safety and performance specifications established by the OEM and to meet its original intended use, any change to the intended use should

<sup>20</sup> 21 CFR 820.3(w).

<sup>21</sup> See the FDA guidance “Reprocessing Medical Devices in Health Care Settings: Validation Methods and Labeling,” available at: <https://www.fda.gov/regulatory-information/search-fda-guidance-documents/reprocessing-medical-devices-health-care-settings-validation-methods-and-labeling>.



***Contains Nonbinding Recommendations******Draft – Not for Implementation***

be evaluated to determine whether the activity is remanufacturing.<sup>22</sup>

2. **Determine whether the activities, individually and cumulatively, significantly change the safety or performance specifications of a finished device** – Under 21 CFR 820.3(w), remanufacturing includes activities that significantly change the performance or safety specifications of the finished device. FDA considers “change” to also include activities that improve the device. Activities that are not *intended to* significantly change the performance or safety specifications, however, should still be evaluated to determine whether they *do* significantly change the finished device’s performance and safety specifications. Multiple changes, when considered cumulatively, may significantly change the performance or safety specifications of the legally marketed device and should be evaluated.
3. **Evaluate whether any changes to a device require a new marketing submission** – Regardless of whether changes made to a legally marketed device are remanufacturing, such changes should be evaluated to determine whether a premarket notification (510(k)) or other marketing submission is required pursuant to the FD&C Act and applicable regulations, and entities should consult relevant guidance for FDA’s recommendations on the topic.<sup>23</sup> For example, a change to a device subject to 510(k) and/or special controls should be considered with respect to the criteria in 21 CFR 807.81 describing when a new 510(k) submission is required and any special controls under the relevant device classification regulation, respectively.
4. **Assess component/part/material<sup>24</sup> dimensional and performance specifications** – Assessment of changes to dimensional and performance specifications can inform whether the activity performed is remanufacturing. The impact of component/part/material changes can be evaluated by comparison to the OEM components/parts/materials specifications and/or through verification and validation testing. Deviations in component/part/material specifications from the OEM’s legally marketed device may result in significant changes to the device’s performance or safety specifications, and may necessitate closer evaluation (such as conducting testing or a risk-based assessment) and consideration of the regulatory criteria describing when a new marketing submission is required.

<sup>22</sup> Consistent with Guiding Principle 3, any changes that affect or change intended use should be considered pursuant to applicable regulations.

<sup>23</sup> See, e.g., “Deciding When to Submit a 510(k) for a Change to an Existing Device,” available at <https://www.fda.gov/regulatory-information/search-fda-guidance-documents/deciding-when-submit-510k-change-existing-device>, and “Deciding When to Submit a 510(k) for a Software Change to an Existing Device,” available at <https://www.fda.gov/regulatory-information/search-fda-guidance-documents/deciding-when-submit-510k-software-change-existing-device>, for FDA’s current thinking on this topic.

<sup>24</sup> 21 CFR 820.3(c) defines a component as any raw material, substance, piece, part, software, firmware, labeling, or assembly which is intended to be included as part of the finished, packaged, and labeled device. In this guidance, “component” and “component/part/material” are used interchangeably. Due to the nature of software and firmware, consideration of whether activities involving them may be remanufacturing is discussed separately from components/parts/materials.

**Contains Nonbinding Recommendations****Draft – Not for Implementation**

- 201 5. **Employ a risk-based approach** – Entities should employ a risk-based approach, such as  
202 one that conforms to or is consistent with ISO 14971: *Medical devices – Application of*  
203 *risk management to medical devices* when assessing whether an activity they perform is  
204 remanufacturing. For the purposes of this guidance, a risk-based assessment is based on  
205 the combination of multiple risk concepts that are important for managing the risks of  
206 medical devices. Risk estimation, risk acceptability, risk control, benefit/risk analysis,  
207 assessment of hazards and hazardous situations, and overall risk evaluation are all  
208 concepts that can be applied during these activities. The concept of risk, as defined in  
209 ISO 14971, is the combination of the probability of occurrence of harm and the severity  
210 of that harm. Although the risk terminology used in this document is primarily derived  
211 from ISO 14971, we recognize that an individual entity’s terminology may differ.  
212

213 For the purposes of this guidance, a new risk is a new hazard or hazardous situation that  
214 did not exist for the legally marketed device. An activity performed on a device may  
215 introduce a new risk, or may modify the probability or severity of a known risk. An  
216 activity is likely remanufacturing when a risk-based assessment identifies any new risks  
217 or significant modifications to known risks, as these are likely to significantly change  
218 performance or safety specifications, in comparison to the legally marketed device.  
219

- 220 6. **Adequately document decision-making** – When deciding whether an activity is  
221 remanufacturing or not, FDA recommends that the rationale for the determination be  
222 documented in sufficient detail, including reference to supporting verification and  
223 validation data, to explain how the determination was made. Specifically, such  
224 documentation should specify why the activities performed on the device do or do not  
225 significantly change the performance or safety specifications, or intended use of the  
226 legally marketed device. If an entity previously determined that an activity was not  
227 remanufacturing, and the same entity is performing the identical activity on the same  
228 version or model of a device, such documentation could reference previous  
229 determinations. Effective documentation can facilitate sound decision-making and  
230 evaluation of relevant factors and information such as adverse events, and provide  
231 important information for an entity to help justify their decision-making in the event that  
232 an inspection is conducted by FDA or this information is otherwise requested.  
233

## 234 **VI. Relevant considerations to determine if activities are** 235 **remanufacturing**

### 236 **A. What is a significant change to device performance or** 237 **safety specifications?**

238 Remanufacturing is the processing, conditioning, renovating, repackaging, restoring, or any other  
239 act done to a finished device that significantly changes the finished device’s performance or  
240 safety specifications, or intended use.<sup>25</sup> For purposes of this draft guidance, FDA generally

---

<sup>25</sup> 21 CFR 820.3(w).

***Contains Nonbinding Recommendations******Draft – Not for Implementation***

considers a significant change to device performance or safety specifications to be one that, based on verification and validation testing and/or a risk-based assessment, results in a finished device that is outside the OEM's performance or safety specifications or introduces new risks or significantly modifies existing risks. For example, a change to a material that contacts the human body and impacts the adequacy of the OEM's validated reprocessing instructions is likely a significant change to device performance or safety specifications, and therefore, is likely remanufacturing. Conversely, replacing an internal capacitor with one that has the same specifications (e.g., same capacitance, working voltage, temperature range, and footprint) is not likely to significantly change device performance or safety specifications and therefore, is likely not remanufacturing.

FDA has identified certain types of activities that, in general, the Agency believes significantly change the legally marketed device's performance or safety specifications:

- Changes to the device's sterilization methods;
- Changes to the device's reprocessing instructions;<sup>26</sup> and
- Changes to the device's control mechanism,<sup>27</sup> operating principle,<sup>28</sup> or energy type.<sup>29</sup>

As discussed below in Section VI.B, activities that result in these changes are likely remanufacturing, and evaluation using the flowchart and accompanying text is not recommended.

Remanufacturing also includes significant changes to a device's intended use (e.g., changing a single-use device to become reusable, changing the anatomical location of use).<sup>30</sup> Therefore, as discussed in Guiding Principle 1, any change to the intended use should be evaluated to determine whether the activity is remanufacturing.

## **B. Determining whether activities are “remanufacturing”**

For activities involving components/parts/materials, FDA recommends the use of the flowchart in this section (Figure 1) to help entities determine if their activities are likely remanufacturing. Although the servicing and remanufacturing definitions and guiding principles in this document apply to software, the approach described in this section should not be applied to software due to

<sup>26</sup> See the FDA guidance “Reprocessing Medical Devices in Health Care Settings: Validation Methods and Labeling,” available at: <https://www.fda.gov/regulatory-information/search-fda-guidance-documents/reprocessing-medical-devices-health-care-settings-validation-methods-and-labeling>.

<sup>27</sup> For purposes of this guidance, a control mechanism is the manner by which the actions of a device are directed. One example of a control mechanism change would be a change from analog to digital control of a medical device.

<sup>28</sup> For purposes of this guidance, an operating principle is the mode of operation or mechanism of action through which a device fulfills (or achieves) its intended use. An example of a new operating principle would be changing the image reconstruction algorithm used in a computed tomography x-ray system from simple back projection to a new, more radiation-efficient method.

<sup>29</sup> For purposes of this guidance, energy type is the type of power input to or output from the device. These changes include both energy output and input changes. A change from emitting microwave energy to radiofrequency (RF) energy would be an example of an energy output change; this type of change would likely be part of a significant redesign.

<sup>30</sup> 21 CFR 820.3(w).

*Contains Nonbinding Recommendations**Draft – Not for Implementation*

its nature and the methods used to evaluate changes to software. Instead, see Section VII for a discussion of changes involving software.

Figure 1 is a visual aid intended to be used in conjunction with the accompanying text in this section and guiding principles. Figure 1 and the accompanying text in this section are intended to address the most common and important considerations that should be evaluated, but is not meant to capture *all* potential considerations that an entity should evaluate to determine if their activities are likely remanufacturing. Rather, they are intended to guide entities in determining when they should further evaluate such activities by conducting testing or a risk-based assessment. Figure 1 and the accompanying text are intended to be consulted after it is determined that there is no significant change to intended use.

In Figure 1, each change (e.g., physical change or change to safety control) should first be assessed individually to determine whether the activity is likely remanufacturing. After evaluating each change individually, the cumulative effects should be assessed to determine whether the activities resulting in the collective changes are likely remanufacturing. The legally marketed device should be used as the basis for comparison for individual changes and the cumulative effects of such changes. When there are no deviations in component/part/material dimensional or performance specifications, or intended use, from the OEM counterpart, and there are no new or modified risks or change in the performance or safety specifications of the legally marketed device, there would likely be no significant changes to the legally marketed device, in the absence of other changes.

FDA does not recommend evaluation with Figure 1 when an activity is performed on behalf of, or otherwise explicitly authorized by, the OEM and the activity returns the legally marketed device to its original performance and safety specifications, and intended use. FDA believes such activities would likely not be remanufacturing, and the determination should be adequately documented.

Entities performing activities on devices should make a determination about whether each activity and the cumulative effects of such activities are remanufacturing and document their rationale.<sup>31</sup> When deciding whether an activity is remanufacturing, entities should document the decision-making process and the basis for the determination. The documentation should be prepared in a way that an FDA investigator or other third party can understand what the change was and the rationale underlying the conclusion. For this, we recommend that the documentation include, at a minimum, the following:<sup>32</sup>

- Product name (including model number and serial number, if applicable);
- Date of activities performed, assessment, and determination;
- Description of device;

<sup>31</sup> In addition, FDA notes that under 21 CFR Part 820, manufacturers are required to maintain certain records as applicable, e.g., service reports.

<sup>32</sup> Consistent with Guiding Principle 6, if the identical activity was previously determined to not be remanufacturing, is being performed by the same entity, and is being performed on the same version or model of a device, such documentation could reference previous determinations.

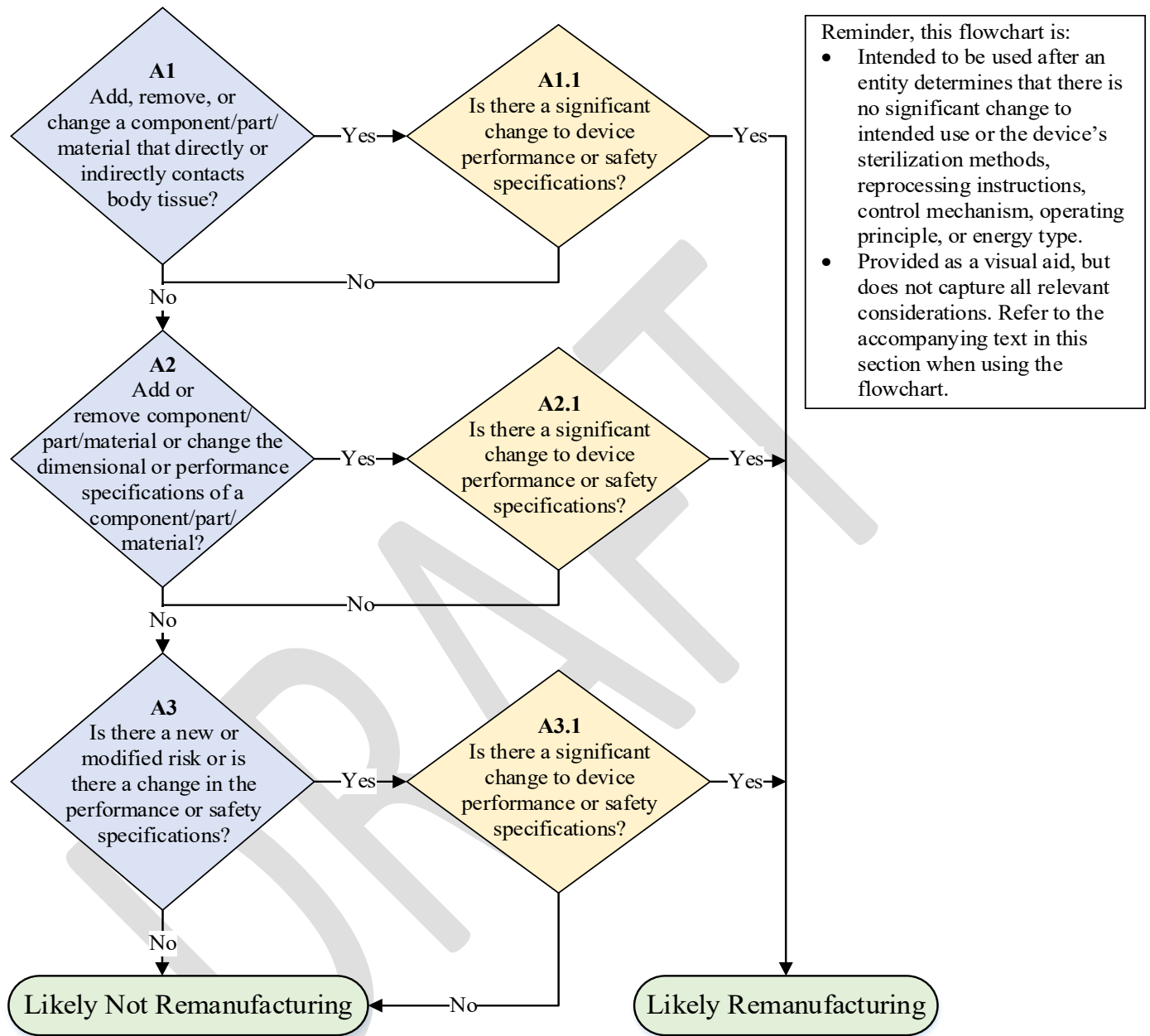
***Contains Nonbinding Recommendations***

***Draft – Not for Implementation***

- Description of activities to be performed, including documentation of components/parts/materials involved;
- Determination of whether the activity is remanufacturing (we recommend using the applicable sections of this guidance);
- Reference to related documents supporting the decision-making process; and
- Signature(s).

FDA has included examples of such documentation in Appendix B.

DRAFT

**Contains Nonbinding Recommendations*****Draft – Not for Implementation***

319 **Figure 1.** Flowchart to help determine whether activities performed are likely remanufacturing.



***Contains Nonbinding Recommendations******Draft – Not for Implementation*****A1. Add, remove, or change a component/part/material that directly or indirectly contacts body tissue?**

Consistent with FDA’s guidance documents on reprocessing<sup>33</sup> and biocompatibility,<sup>34</sup> respectively, entities should assess how their activities may affect validated reprocessing instructions or cause an unacceptable adverse biological response resulting from device contact with the human body, including both patient and healthcare provider tissue.

Direct contact is when a component/part/material comes into physical contact with body tissue, such as catheter tubing used on a patient. A component/part/material has indirect contact when a fluid or gas passes through it prior to the fluid or gas coming into physical contact with body tissue (i.e., the device or component/part/material itself does not physically contact body tissue). For example, materials in a catheter hub (the part of the catheter that is external to the patient) indirectly contact the patient when fluids or drugs are infused through the hub and into the patient. Both direct and indirect contact with the patient or user of the device should be considered when answering A1.

If there is any addition, removal, or change to a component/part/material on the finished device, and that component/part/material directly or indirectly contacts body tissue, the answer to A1 should be “yes.” This includes exposing a previously unexposed component/part/material to direct or indirect contact with body tissue. Additionally, if there is any change in material type, formulation, or chemical composition for a component/part/material that directly or indirectly contacts body tissue, the answer to A1 should be “yes.” If the entity is uncertain how to respond to A1, the answer should be “yes.” A “yes” answer to A1 does not necessarily mean that the activity is remanufacturing. Rather, when an entity makes such changes, it should analyze the impact of the change on the device’s performance and safety specifications using the text in A1.1.

If no component/part/material added, removed, or changed directly or indirectly contacts body tissue, the answer should be “no” and then proceed to A2.

**A1.1 Is there a significant change to device performance or safety specifications?**

If the activity adds, removes, or changes a component/part/material that directly or indirectly contacts body tissue (as mentioned above, this includes an activity that exposes a previously unexposed component/part/material to body tissue either directly or indirectly), a risk-based assessment should be conducted. The assessment should be conducted to determine whether there is a significant change to the biocompatibility or the validated reprocessing instructions of

<sup>33</sup> See the FDA guidance “Reprocessing Medical Devices in Health Care Settings: Validation Methods and Labeling,” available at: <https://www.fda.gov/regulatory-information/search-fda-guidance-documents/reprocessing-medical-devices-health-care-settings-validation-methods-and-labeling>.

<sup>34</sup> See the FDA guidance “Use of International Standard ISO 10993-1, ‘Biological evaluation of medical devices – Part 1: Evaluation and testing within a risk management process,’” available at: <https://www.fda.gov/regulatory-information/search-fda-guidance-documents/use-international-standard-iso-10993-1-biological-evaluation-medical-devices-part-1-evaluation-and>.

***Contains Nonbinding Recommendations******Draft – Not for Implementation***

the legally marketed device. An activity that results in such change may be considered remanufacturing.

Depending on the magnitude of the change and the nature of the component/part/material, reprocessing validation and a comprehensive biocompatibility risk assessment or testing may be necessary. Entities should incorporate factors that affect the reprocessing and biocompatibility of a device in their risk-based assessment and testing where appropriate. These factors may include the materials of construction, the processing of the materials, methods (including the sterilization process), any residuals from aids used during the process, and intended use life of the legally marketed device. Activities that impact the adequacy of the legally marketed device's validated reprocessing instructions are likely remanufacturing.

If the answer to A1.1 is "yes," then the activity would likely be remanufacturing. If the answer to A1.1 is "no," then proceed to A2.

## **A2. Add or remove component/part/material or change the dimensional or performance specifications of a component/part/material?**

*Add or remove component/part/material?* If there is any addition of a component/part/material to a legally marketed device that was not originally part of the legally marketed device, the answer to A2 should be "yes." Examples include adding an adhesive to mend a break in the device or fasteners to secure a component/part/material. If there is any removal of a component/part/material to a legally marketed device that is not replaced in the legally marketed device, the answer to A2 should be "yes." Examples include removing a fastener or barrier without replacement. Add or remove component/part/material also includes replacing an OEM component/part/material with the same OEM component/part/material or a non-OEM component/part/material.<sup>35</sup>

*Change the dimensional or performance specifications of a component/part/material?* If there is any change to or replacement of a component/part/material of the legally marketed device, which affects the component/part/material's dimensional or performance specifications, the answer to A2 should be "yes."

If a component/part/material is not being added or removed, or the dimensional or performance specifications of a component/part/material are not being changed, the answer to A2 should be "no." If uncertain, the answer to A2 should be "yes."

A "yes" answer to A2 does not necessarily mean that the activity is remanufacturing. Rather, when an entity makes such changes, it should analyze the impact of the change on the device's performance and safety specifications using the text in A2.1. If the answer to A2 is "no," then proceed to A3.

<sup>35</sup> As discussed above in Section VI.B., FDA does not recommend evaluation with Figure 1 when an activity is performed on behalf of, or otherwise explicitly authorized by, the OEM and the activity returns the legally marketed device to its original performance and safety specifications, and intended use. FDA believes such activities would likely not be remanufacturing, and the determination should be adequately documented.

**Contains Nonbinding Recommendations****Draft – Not for Implementation****A2.1 Is there a significant change to device performance or safety specifications?**

*Does the added or removed component/part/material significantly change the device performance or safety specifications?* When evaluating whether the addition or removal of a component/part/material significantly changes the device's performance or safety specifications, the entity should consider the intended use life of the legally marketed device. For instance, many reusable devices are reprocessed numerous times within their intended use life. Applicable considerations should include an assessment of whether the added component will withstand repeated reprocessing cycles within the device's intended use life or whether the removed component exposes previously unexposed components that will withstand repeated reprocessing cycles within the device's intended use life. Such an assessment can include verification and validation testing or a risk-based assessment describing why such testing is not warranted. If the reusable device will not be able to withstand repeated reprocessing cycles within its intended use life, the addition or removal of the component/part/material may significantly change the legally marketed device's performance or safety specifications.

*Do the changed dimensional specifications of the component/part/material significantly change the device performance or safety specifications?* In determining whether an activity is remanufacturing for these types of changes, the entity should consider not only the magnitude of the dimensional specification change, but the criticality of the modified dimension. The entity should consider whether dimensional specifications meet a minimum or maximum specification (e.g., outer diameter cannot exceed 3.0 mm) or are within a range of acceptable tolerance specifications. If dimensional specifications are within the acceptable range, the answer to A2.1 would likely be "no;" however, for changes that are outside the acceptable range of dimensional specifications, the answer to A2.1 would likely be "yes."

*Do the changed performance specifications of the component/part/material significantly change the device performance or safety specifications?* When evaluating if there is a significant change to performance or safety specifications, the entity should consider whether performance outputs meet a minimum and/or maximum specification (e.g., temperature within chamber cannot exceed 25 °C and pressure cannot be less than 150 kPa) or are within a range of acceptable tolerance specifications (e.g., pump flowrate must be between 2 and 20 mL/hour). If performance specifications are within the acceptable range, the answer to A2.1 would likely be "no;" however, for changes that result in performance specifications that are outside the acceptable range, the answer to A2.1 would likely be "yes."

If the answer to A2.1 is "yes," then the change would likely be remanufacturing. If the answer to A2.1 is "no," then proceed to A3.

**A3. Is there a new or modified risk or is there a change in the performance or safety specifications?**

The entity should perform a risk-based assessment to identify new or modified risks or a change in the performance or safety specifications of the legally marketed device based on the activity being performed on the device. Both the individual change and cumulative changes performed

***Contains Nonbinding Recommendations***

***Draft – Not for Implementation***

on the legally marketed device should be considered. While individual changes may not significantly change the legally marketed device’s performance or safety specifications, the cumulative changes may do so. The extent of the assessment should be appropriate considering the nature and extent of the activities being performed.

*Is there a new or modified risk?* A risk-based assessment can identify whether there are new risks or modified existing risks in comparison to the legally marketed device. If a new risk is created or an existing risk has been modified based on the activity being performed, the answer to A3 should be “yes,” and this activity should be evaluated using the text in A3.1. If uncertain, the answer to A3 should be “yes.”

*Is there a change in the performance or safety specifications?* A risk-based assessment can also identify whether there is a change in performance or safety specifications. This assessment should consider, for example, how a change could impact a device’s continued conformity to a voluntary consensus standard or compliance with a regulation, such as special controls identified in a device classification regulation. This assessment should also consider whether activities that break a seal or barrier can adequately return the device to its legally marketed performance and safety specifications. If a change to performance or safety specifications has been identified, the answer to A3 should be “yes.” If uncertain, the answer to A3 should be “yes.”

When an entity makes a change that has a “yes” answer to A3, the entity should analyze the impact of the change on the device’s performance and safety specifications using the text in A3.1. If the answer to A3 is “no,” then the change is likely not remanufacturing.

### **A3.1 Is there a significant change to device performance or safety specifications?**

If new or modified risks were identified, the entity should evaluate whether they significantly change the legally marketed device’s performance or safety specifications using the output of the risk-based assessment performed in A3. Altering or bypassing a safety feature (e.g., fuses, alerts, alarms, interlocks) likely significantly changes the legally marketed device’s performance or safety specifications. Changes that impact compliance with a regulation or alter conformity with a voluntary consensus standard would likely significantly change the legally marketed device’s performance or safety specifications and may also adulterate and/or misbrand the device.<sup>36</sup>

If the answer to A3.1 is “yes,” then the change would likely be remanufacturing. If the answer to A3.1 is “no,” then the change is likely not remanufacturing.

## **VII. Changes involving software**

As described in Section VI, Figure 1 and its accompanying text should not be applied to changes involving software. Many software changes are likely remanufacturing because of their impact on a product’s software architecture, software requirements specifications, unresolved anomalies,

<sup>36</sup> See, e.g., sections 501(e)(2) and 502(o) of the FD&C Act.

***Contains Nonbinding Recommendations******Draft – Not for Implementation***

and other key characteristics. Further, because the probability of a software failure cannot be determined using traditional statistical methods, the risk-based assessment approach that FDA recommends in Section VI should not be applied to software changes. Instead, FDA has identified certain activities performed on software that are likely not remanufacturing because they generally do not significantly change the performance or safety specifications of the device:

- Activities performed on behalf of or otherwise explicitly authorized by the OEM that return the legally marketed device to its performance and safety specifications, and intended use;
- Implementing OEM provided updates and upgrades;
- Running software-based hardware diagnostics;
- Assessing for viruses, malware, and other cybersecurity related issues;
- Reinstalling OEM software to restore original performance and safety specifications;
- Reverting software to a previous configuration;
- Installing cybersecurity updates that are authorized by the OEM;
- Turning on or off connectivity features (e.g., Wi-Fi and Bluetooth connections) consistent with OEM intended use;
- Performing data backup and recovery operations;
- Assessing software inventory;
- Collecting system logs;
- Managing user accounts; and
- Accessing diagnostic and repair information.

Other activities involving changes to software are likely to significantly change a device's performance or safety specifications, such that the activity is likely remanufacturing. However, if an entity believes that an activity involving a change to software does not significantly change a device's performance or safety specifications, the entity should adequately document its decision-making (see Guiding Principle 6). Any activity involving software changes that significantly modifies a device's intended use would be remanufacturing.<sup>37</sup>

Entities should also consider the unintended consequences and cumulative effects of any software change(s). Entities performing activities on devices should make a determination about whether each activity and the cumulative effects of the changes resulting from the activities are remanufacturing and document their rationale.

## **VIII. Considerations for labeling**

Based on publicly available information and FDA's activities discussed in Section II.A of this draft guidance, FDA believes that OEMs of reusable devices intend for their devices to routinely undergo both preventive maintenance and repair. It is important that such devices include instructions on how to adequately return a device to its performance and safety specifications

---

<sup>37</sup> See 21 CFR 820.3(w).



***Contains Nonbinding Recommendations******Draft – Not for Implementation***

established by the OEM.<sup>38</sup> Unintentional remanufacturing can occur when entities do not have the instructions necessary to return a device to its original performance and safety specifications. The lack of adequate servicing instructions can also create challenges in the availability of quality, safe, and effective devices.

Consistent with promoting and protecting the public health, FDA encourages OEMs, as an industry best practice, to provide servicing instructions that facilitate routine maintenance and repair of their reusable devices.<sup>39</sup> FDA recommends that the labeling of reusable devices include the following information, as applicable, to facilitate routine device maintenance and repair:

- A description of the key performance and safety specifications;
- Critical technical or functional specifications, including:
  - Physical dimensions;
  - Electrical characteristics, including batteries (e.g., chemistry, amperage, voltage, rechargeability), internal fuses, and power supply (e.g., voltage, amperage, frequency); and
  - Device-specific performance specifications (e.g., flow rate accuracy or range, humidity, temperature, wavelength).
- The recommended maintenance activities and schedule;
- Recommended routine testing and acceptance criteria to confirm that the device remains within its performance and safety specifications;
- A description of error codes, alerts, and alarm features on the device;
- Precautions and warnings relevant to servicing the device; and
- Version number and release date of software.

---

<sup>38</sup> Section 502(f)(1) of the FD&C Act requires that labeling bear adequate directions for use. For non-prescription devices, adequate directions for use include instructions on preparing a device for use. 21 CFR 801.5(g). Prescription devices are exempt from the adequate directions for use requirement provided certain conditions are met, including that the labeling bear “information for use, including indications, effects, routes, methods, and frequency and duration of administration, and any relevant hazards, contraindications, side effects, and precautions under which practitioners licensed by law to administer the device can use the device safely and for the purpose for which it is intended...” 21 CFR 801.109(c).

<sup>39</sup> FDA’s recommendations in this guidance are not intended to encourage the disclosure of trade secrets or confidential commercial information (CCI).



**Contains Nonbinding Recommendations****Draft – Not for Implementation****Appendix A. Examples**

The following are illustrative examples of activities that may be performed on devices with explanations about why such examples are or are not likely remanufacturing. Note that these generalized examples do not necessarily account for every possible detail, risk, or consideration that a manufacturer should evaluate, and should not be taken to mean that the changes described are or are not definitively remanufacturing. Real-world decisions will depend on the specific facts and circumstances, including the specific details of the changes made to the specific device.

**(1) Component/part/material activities****Example E.1**

**Activity:** The door of an infusion pump was bent and now pinches the administration set. The flow rate accuracy fell outside the OEM's specified accuracy range. The door is replaced with a non-OEM door that is marketed as compatible with this infusion pump. It has the same overall dimensions and is made from a similar material of construction. However, the replacement door material is more rigid than the original door.

**Relevant questions:**

*A1. Add, remove, or change a component/part/material that directly or indirectly contacts body tissue?*

No. The existing and replacement doors do not have direct or indirect contact with the patient's body tissue.

*A2. Add or remove component/part/material or change the dimensional or performance specifications of a component/part/material?*

Yes, the old door was removed and replaced. While the new door is marketed as compatible, all dimensions were confirmed through comparative measurement, including the hinges and latch. The specific material of the original door is unknown and there is a noticeable difference in flexibility that may impact the pump's performance specifications.

*A2.1 Is there a significant change to device performance or safety specifications?*

No. Once replaced, the door was confirmed to open and close with similar effort as the original door and it was confirmed that the added rigidity did not significantly change the pump's performance or safety specifications (e.g., flowrate accuracy).

*A3. Is there a new or modified risk or is there a change in the performance or safety specifications?*

No. A risk-based assessment determined that there are no new or modified risks and there is no change in performance or safety specifications (e.g., the change does not alter conformity to a voluntary consensus standard or compliance with a regulation).

**Decision:** Not Remanufacturing.

***Contains Nonbinding Recommendations******Draft – Not for Implementation*****Example E.2**

**Activity:** The rotor within a peristaltic infusion pump no longer functions as intended and is replaced. The subject pump rotor is no longer supported by the OEM, but a comparable off-the-shelf rotor is available. The dimensions of the rotor, including the individual rollers, are the same; however, the material of construction of the rollers, which contact and apply pressure to the administration set, appears to be stainless steel. This is different from the plastic rollers in the legally marketed device.

**Relevant questions:**

*A1. Add, remove, or change a component/part/material that directly or indirectly contacts body tissue?*

No. Neither the existing or replacement component directly or indirectly contact body tissue. It is only in contact with the outside of the administration set.

*A2. Add or remove component/part/material or change the dimensional or performance specifications of a component/part/material?*

Yes. The rotor was removed and replaced. Also, although the dimensional specifications of the non-OEM pump rotor, including the individual rollers, are the same as the OEM rotor, the roller materials are different.

*A2.1 Is there a significant change to device performance or safety specifications?*

Yes. Once the rotor was replaced, the device appears to function adequately. The change in material of the rollers does not significantly change the accuracy of the flowrate across the labeled flowrate range. However, a risk-based assessment identified that the change in material of the rollers can affect the useful life of the administration set. The change in the roller material from plastic to stainless steel increases the administration set wear and/or breakage due to fatigue. Evaluation of this risk concluded that the increased fatigue on the administration set is more likely to lead to patient under-dosing before the administration set is intended to be replaced. This significantly changes the device's performance and safety specifications.

**Decision:** Remanufacturing.

**Example E.3**

**a. Activity:** The gradient coil of a magnetic resonance (MR) system was damaged during an imaging session and needs to be replaced. The gradient coil is replaced with a non-OEM gradient coil. The maximum slew rate of the coil matches that of the OEM gradient coil; however, the peak gradient strength is larger than the OEM coil.

**Relevant questions:**

*A1. Add, remove, or change a component/part/material that directly or indirectly contacts body tissue?*

No. The gradient coil does not have direct or indirect contact with body tissue.

## ***Contains Nonbinding Recommendations***

### ***Draft – Not for Implementation***

*A2. Add or remove component/part/material or change the dimensional or performance specifications of a component/part/material?*

Yes. The gradient coil was removed and replaced, and the new gradient coil has a larger peak gradient strength.

*A2.1 Is there a significant change to device performance or safety specifications?*

Yes. An assessment was performed to determine the significance of the change. A gradient coil with a larger peak gradient strength significantly changes the imaging performance specifications (e.g., slice thickness, spatial resolution).

**Decision:** Remanufacturing.

- b. Activity:** The gradient coil of an MR system was damaged during an imaging session and needs to be replaced. It is replaced with a non-OEM gradient coil that has different dimensional specifications and coil design.

**Relevant questions:**

*In this example, the answers to flowchart questions A1 and A2 are the same as Example E.3.a. except that for A2, the new gradient coil has different dimensional specifications and coil design.*

*A2.1 Is there a significant change to device performance or safety specifications?*

No. The new gradient coil only differs by small changes in design and dimensional specifications. There are no significant changes to the performance and safety specifications (e.g., slew rate, peak gradient strength, power).

*A3. Is there a new or modified risk or is there a change in the performance or safety specifications?*

No. A risk-based assessment identified no new or modified risks or change in the performance or safety specifications due to this change because the non-OEM gradient coil has the same hardware performance specifications (e.g., slew rate), equivalent imaging performance, and meets the same safety and performance specifications (e.g., acoustic output) when compared to the OEM gradient coil.

**Decision:** Not Remanufacturing.

#### **Example E.4**

**Activity:** The slide heater pads on an immunohistochemistry (IHC) autostainer are worn out and need to be replaced. They are replaced with an OEM part.

**Relevant questions:**

*A1. Add, remove, or change a component/part/material that directly or indirectly contacts body tissue?*

No. The slide heater pads do not have direct or indirect contact with body tissue.

## ***Contains Nonbinding Recommendations***

### ***Draft – Not for Implementation***

*A2. Add or remove component/part/material or change the dimensional or performance specifications of a component/part/material?*

Yes. The heater pad components were physically removed and replaced with new pads.

*A2.1 Is there a significant change to device performance or safety specifications?*

No. An assessment was performed to evaluate this replacement and identified no changes to dimensions, materials, or performance or safety specifications of the pads.

*A3. Is there a new or modified risk or is there a change in the performance or safety specifications?*

No. A risk-based assessment identified no new or modified risks because the slide heater pads are identical to the original part from the OEM. The device now functions within its functional specifications identified in the labeling. There is no change in the performance or safety specifications.

**Decision:** Not Remanufacturing.

#### **Example E.5**

**Activity:** The tubing on a sample processor became kinked from use and needs to be replaced. Tubing was found from the same OEM but the tubing is intended for use with a different sample processor.

#### **Relevant questions:**

*A1. Add, remove, or change a component/part/material that directly or indirectly contacts body tissue?*

No. There is no direct or indirect contact between the tubing and body tissue.

*A2. Add or remove component/part/material or change the dimensional or performance specifications of a component/part/material?*

Yes. The tubing was removed and replaced with new tubing of a different inner diameter.

*A2.1 Is there a significant change to device performance or safety specifications?*

Yes. The inner diameter of the tubing is different from the legally marketed device. Verification and validation testing was performed to evaluate this replacement and identified significant changes to performance because different fluid characteristics (e.g., flow rate) than those specified for the legally marketed device were noted with the new tubing.

**Decision:** Remanufacturing.

#### **Example E.6**

**Activity:** A tissue pre-treatment water bath was updated by replacing the heating chamber with one that has a different temperature range.

## ***Contains Nonbinding Recommendations***

### ***Draft – Not for Implementation***

#### **Relevant questions:**

*A1. Add, remove, or change a component/part/material that directly or indirectly contacts body tissue?*

No. The tissue specimens have been removed from the human body, are within a sealed container, and neither the water bath nor heating chamber directly or indirectly contacts the tissue.

*A2 Add or remove component/part/material or change the dimensional or performance specifications of a component/part/material?*

Yes. The heating chamber was removed and replaced. The heating chamber's performance specifications were changed because the new heating chamber has a different temperature range.

*A2.1 Is there a significant change to device performance or safety specifications?*

Yes. The performance is significantly changed because the heating range extends beyond that of the heating chamber in the legally marketed device.

**Decision:** Remanufacturing.

#### **Example E.7**

- a. **Activity:** A stainless steel manual drill is intended to be used in the implantation of orthopedic devices. The drill is intended to be reprocessed and reused for multiple procedures. The drill was sharpened because it is dull and difficult to use. This is the first time the drill has been sharpened.

#### **Relevant questions:**

*A1. Add, remove, or change a component/part/material that directly or indirectly contacts body tissue?*

Yes. Sharpening the drill removes material and exposes a fresh surface that directly contacts bone.

*A1.1 Is there a significant change to device performance or safety specifications?*

No. The drill is not coated. The material and structure of the drill that contacts body tissue is uniform. A risk-based assessment concluded that removal of material due to sharpening does not significantly change the biocompatibility or reprocessing.

*A2. Add or remove component/part/material or change the dimensional or performance specifications of a component/part/material?*

Yes. Sharpening of the drill removes material changing the dimensions of the drill.

*A2.1 Is there a significant change to device performance or safety specifications?*

No. The drill was returned to its performance and safety specifications because the entity sharpened the device to its labeled outer diameter and original edge profile angle.

***Contains Nonbinding Recommendations******Draft – Not for Implementation***

758 *A3. Is there a new or modified risk or is there a change in the performance or safety*  
759 *specifications?*

760 Yes. Sharpening the drill may change the size of the resulting pilot drill hole. Changing  
761 the size of the pilot hole can change the fit of the implant or overall purchase in bone  
762 such that the mechanical integrity of the implant is compromised.

763  
764 *A3.1 Is there a significant change to device performance or safety specifications?*

765 No. Based on the facility's maintenance record, it was determined that this is the first  
766 drill sharpening. The drill produces the same pilot hole size as the legally marketed  
767 device after the sharpening has been completed. There is no significant change to the  
768 device's performance or safety specifications at this time.

769  
770 **Decision:** Not Remanufacturing.

- 771  
772 b. **Activity:** A stainless steel manual drill with a titanium nitride coating is intended to be  
773 used in the implantation of orthopedic devices. The drill is intended to be reprocessed and  
774 reused for multiple procedures. The drill was sharpened because it is dull and difficult to  
775 use. The drill has been sharpened multiple times.

776  
777 **Relevant questions:**

778 *A1. Add, remove, or change a component/part/material that directly or indirectly*  
779 *contacts body tissue?*

780 Yes. Sharpening the drill removes material and exposes a fresh surface that directly  
781 contacts bone.

782  
783 *A1.1 Is there a significant change to device performance or safety specifications?*

784 No. While sharpening the drill exposes the stainless steel surface beneath the coating,  
785 both the surface coating and underlying stainless steel have been subjected to a  
786 biocompatibility assessment. Additionally, a risk-based assessment concluded that  
787 removal of material due to sharpening does not significantly change the biocompatibility  
788 or reprocessing.

789  
790 *A2. Add or remove component/part/material or change the dimensional or performance*  
791 *specifications of a component/part/material?*

792 Yes. Sharpening of the drill removes material changing the dimensions and cutting  
793 surface of the drill.

794  
795 *A2.1 Is there a significant change to device performance or safety specifications?*

796 Yes. Based on the facility's maintenance record, it was determined that the drill has been  
797 sharpened multiple times. While the outer diameter of the drill is not significantly  
798 changed from the legally marketed device, the titanium nitride coating is no longer intact  
799 on the cutting surface of the drill, causing inefficient or destructive cutting. This activity  
800 significantly changes the device's performance and safety specifications.

801  
802 **Decision:** Remanufacturing.



**Contains Nonbinding Recommendations****Draft – Not for Implementation****Example E.8**

- a. **Activity:** The lens of an endoscope is cracked. The lens is affixed by an epoxy that is not described in the labeling. The cracked lens was removed and replaced. The epoxy used was purchased from the OEM and is identical to that used in the legally marketed device. The replacement lens was not purchased from the OEM. The lens was tested and demonstrated to have the same optical specifications (e.g., focal length, Abbe number) and materials as the original lens.

**Relevant questions:**

*A1. Add, remove, or change a component/part/material that directly or indirectly contacts body tissue?*

Yes, both the lens and the epoxy directly contact body tissue.

*A1.1 Is there a significant change to device performance or safety specifications?*

No. The epoxy is identical to the epoxy used in the legally marketed device. The replacement lens is the same material as original lens. A risk-based assessment that considered both the individual and cumulative changes was performed to determine if the procedure used to replace the lens affects biocompatibility and reprocessing instructions. A biocompatibility assessment confirmed that there are no new surfaces previously unexposed to body tissue. A comprehensive reprocessing risk assessment and testing demonstrated that the validated reprocessing instructions identified in the labeling of the legally marketed device are not impacted by the replacement parts or the procedure used to replace the parts.

*A2. Add or remove component/part/material or change the dimensional or performance specifications of a component/part/material?*

Yes. The epoxy and lens were replaced.

*A2.1 Is there a significant change to device performance or safety specifications?*

No. The optical performance testing (e.g., resolution and distortion) and reprocessing risk assessment and testing indicated there has been no significant change in performance or safety specifications.

*A3. Is there a new or modified risk or is there a change in the performance or safety specifications?*

No. A risk-based assessment was performed that considered both the individual and cumulative changes that could have affected biocompatibility, reprocessing, and optical performance. This assessment identified that there are no new or modified risks, and there is no change in performance or safety specifications.

**Decision:** Not Remanufacturing.

- b. **Activity:** The lens of an endoscope is cracked. The lens is affixed by an epoxy that is not described in the labeling. The cracked lens was removed and replaced. The epoxy used

**Contains Nonbinding Recommendations*****Draft – Not for Implementation***

was purchased from the OEM and is identical to that used in the legally marketed device. The replacement lens comes from a different endoscope model from the same OEM; that model was 510(k)-cleared with improved optical performance (e.g., resolution and distortion) relative to the original endoscope. The replacement lens has the same material but different optical specifications (e.g., focal length, Abbe number) from the original.

**Relevant questions:**

*In this example, the answers to flowchart questions A1, A1.1, and A2 are the same as Example E.8.a.*

*A2.1 Is there a significant change to device performance or safety specifications?*

Yes. The epoxy is identical to that used in the legally marketed device, but the lens has different optical specifications from the original lens. The endoscope with the replacement lens has different imaging specifications relative to the legally marketed device. While the replacement lens is present on another 510(k)-cleared device, it was not present on the original endoscope and significantly changes the performance specifications of the original endoscope.

**Decision:** Remanufacturing.

**Example E.9**

**Activity:** An endoscope's connection to the video processor was damaged during use. After repair, it was observed that the endoscope readily disconnected from the video processor. To address this problem, an adapter was added to reduce the probability of a disconnection between the endoscope and video processor. The adapter was found to be capable of connecting to the video processor; however, it is bulkier than the connector.

**Relevant questions:**

*A1. Add, remove, or change a component/part/material that directly or indirectly contacts body tissue?*

No. The added adapter does not directly or indirectly contact body tissue.

*A2. Add or remove component/part/material or change the dimensional or performance specifications of a component/part/material?*

Yes, the adapter has been added to the endoscope.

*A2.1 Is there a significant change to device performance or safety specifications?*

No. The adapter still allows the endoscope to be connected to the video processor and optical performance testing demonstrated the same optical performance as the original endoscope.

*A3. Is there a new or modified risk or is there a change in the performance or safety specifications?*

Yes. A risk-based assessment was performed to determine the effects of this added component. Increased risks exist with the added adapter, such as disconnection from the light

***Contains Nonbinding Recommendations******Draft – Not for Implementation***

source, and the potential change to the electrical safety and electromagnetic compatibility (EMC) of the device.

*A3.1 Is there a significant change to device performance or safety specifications?*

Yes. Disconnection from a light source during a procedure could result in a loss of imaging and adverse events such as increased procedure time or other patient injuries such as perforation. Additionally, testing should also be performed for the electrical safety and EMC of the device.

**Decision:** Remanufacturing.

**Example E.10**

**Activity:** The motor on a powered wheelchair no longer functions and does not propel the wheelchair as intended. The motor was inspected and it was determined that the motor should be replaced. Neither the identical motor nor one with similar specifications could be located. A motor of similar size was inserted with different power and speed specifications.

**Relevant questions:**

*A1. Add, remove, or change a component/part/material that directly or indirectly contacts body tissue?*

No. The motor does not directly or indirectly contact body tissue.

*A2. Add or remove component/part/material or change the dimensional or performance specifications of a component/part/material?*

Yes. The original motor was removed and replaced.

*A2.1 Is there a significant change to device performance or safety specifications?*

Yes. While the motor has the same physical dimensions, the replacement motor has a different power output and maximum speed than the legally marketed device. This significantly changes the device's performance specifications because the wheelchair can go faster than intended. This also significantly changes the device's safety specifications because the controller and software to operate the wheelchair may no longer be compatible with the motor.

**Decision:** Remanufacturing.

**Example E.11**

- a. **Activity:** The liquid cooling system responsible for maintaining the temperature of a transcranial magnetic stimulation (TMS) coil is malfunctioning and causing the system to overheat. The cooling system was inspected and it was determined that the pump circulating the liquid coolant stopped functioning. A replacement pump was located and installed with no additional changes to the device.

**Contains Nonbinding Recommendations****Draft – Not for Implementation****Relevant questions:**

*A1. Add, remove, or change a component/part/material that directly or indirectly contacts body tissue?*

No. The liquid coolant is maintained in the sealed coolant system and neither the liquid coolant nor the pump directly or indirectly contacts body tissue.

*A2. Add or remove component/part/material or change the dimensional or performance specifications of a component/part/material?*

Yes, the pump was replaced.

*A2.1 Is there a significant change to device performance or safety specifications?*

No. Both the dimensions and performance specifications of the original pump were assessed in comparison to the replacement part. The replacement pump has the same dimensional and performance specifications of the original pump. The overall performance and safety specifications of the TMS coils were verified by testing to be the same.

*A3. Is there a new or modified risk or is there a change in the performance or safety specifications?*

No. A risk-based assessment identified no new or modified risks because the replacement pump is equivalent to that used in the OEM's legally marketed device and there is no change in the device performance or safety specifications.

**Decision:** Not Remanufacturing.

- b. **Activity:** The liquid cooling system responsible for maintaining the temperature of a TMS coil is malfunctioning and causing the system to overheat. The cooling system was inspected and it was determined that the pump circulating the liquid coolant stopped functioning. A replacement pump was located with the same size and flow specifications, but it uses a different coolant liquid. The pump was replaced with one that uses a different coolant into the cooling system.

**Relevant questions:**

*In this example, the answer to flowchart question A1 is the same as Example E.11.a.*

*A2. Add or remove component/part/material or change the dimensional or performance specifications of a component/part/material?*

Yes. A replacement pump that uses a different coolant liquid was installed.

*A2.1 Is there a significant change to device performance or safety specifications?*

Yes. Although the pump has the same dimensional and flow specifications as the original pump, the new pump uses a different liquid coolant. The new liquid coolant does not have the same heat capacity as that used in the legally marketed device. Verification and validation testing was performed and it was determined that there was a significant

**Contains Nonbinding Recommendations**

***Draft – Not for Implementation***

change to cooling effectiveness, which poses a safety hazard when the TMS coil is not properly cooled. This may burn the patient or cause further device malfunctions.

**Decision:** Remanufacturing.

**Example E.12**

**Activity:** An energy-delivering aesthetic device has multiple compatible handpieces with specific areas of application. Applicator A can only be used for the chin, while Applicator B can only be used on the abdomen. An entity cannibalizes Applicator B and uses those parts to repair Applicator A for use on the chin.

**Relevant questions:**

*A1. Add, remove, or change a component/part/material that directly or indirectly contacts body tissue?*

Yes. The distal end of Applicator B is used to reconstruct Applicator A. It directly contacts the patient and delivers the energy.

*A1.1 Is there a significant change to device performance or safety specifications?*

No. The distal end of both applicators has identical materials and the reprocessing instructions provided by the OEM are the same for both applicators. A risk-based assessment was performed to determine the effects of implementing these repairs on the biocompatibility and reprocessing. A biocompatibility assessment and reprocessing risk assessment were used to determine that the performance and safety specifications of the device were not significantly changed.

*A2. Add or remove component/part/material or change the dimensional or performance specifications of a component/part/material?*

Yes. The distal end of Applicator B has different dimensional specifications compared to Applicator A.

*A2.1. Is there a significant change to device performance or safety specifications?*

Yes. The surface area that contacts the patient has increased by 150%. The increase in surface area changes the energy output delivered to the patient, which significantly changes both the performance and safety specifications of Applicator A.

**Decision:** Remanufacturing.

**(2) Software activities**

**Example S.1**

**Activity:** A specular microscope with a camera is intended for examination of corneal endothelium and for measurement of the thickness of the cornea. The software was updated to implement an OEM-authorized patch.

*Contains Nonbinding Recommendations**Draft – Not for Implementation*

**Relevant analysis:** The installation of this OEM-authorized patch does not significantly change the device performance or safety specifications. See Section VII of this draft guidance for further discussion of changes involving software. The patch is intended to maintain the original specifications.

**Decision:** Not Remanufacturing.

**Example S.2**

**Activity:** A device has the capability of real-time remote customer service where the current status of the device can be accessed. A capability is added so that the customer service technician can access and directly manipulate the device, including changing device settings, resetting the device, delivering energy, and positioning the device.

**Relevant analysis:** The capability of the customer service technician to control the device introduces new risks (e.g., accidental device reset, unintended device movement) and functionality (remote control and access) that significantly changes the finished device's performance and safety specifications.

**Decision:** Remanufacturing.

**Example S.3**

**Activity:** A device that connects to a facility's network has software that was designed to run the Microsoft Windows operating system (OS). Adjustments are made to allow the device to run using a Linux OS.

**Relevant analysis:** This change introduces new risks and may impact mitigations for existing risks that significantly change the finished device's performance and safety specifications. This is a redesign of the product and includes the addition of integration with both device drivers for the target OS as well as specific features of the OS.

**Decision:** Remanufacturing.



*Contains Nonbinding Recommendations**Draft – Not for Implementation***Appendix B. Documentation examples**

The examples below are to illustrate one possible approach to documentation; other approaches may also be appropriate. Entities are encouraged to use an approach that works for their specific purposes, taking into account the considerations discussed above. The first example demonstrates a simple change that does not necessitate detailed analysis. The second example demonstrates a more complex change for which additional analysis and reference to supporting documentation are warranted. These are generalized examples to demonstrate documentation principles and do not necessarily account for every possible detail, risk, or consideration.

**Remanufacturing Assessment  
(Example 1)**

**Product:** Pump ABC, Serial# 123-456

**Date of activities performed:** 12/11/2018

**Date assessment performed:** 12/10/2018

**Description of device:** Syringe pump

**Description of activities performed:** Replaced broken door with part #xxx

**Determination of whether the activity is remanufacturing:** While a change to a body contacting component, the door used was OEM-provided and is identical to the broken door. Because it is a replacement of an identical part, there are no changes to performance or safety specifications. This activity is not remanufacturing.

**Reference to related documents supporting the decision-making process:** N/A

**Technician performing service:** xxx

**Reviewed by:** xxx

**Signature(s):** xxx

***Contains Nonbinding Recommendations******Draft – Not for Implementation*****Remanufacturing Assessment  
(Example 2)****Product:** Endoscope Infinity, Serial #4FR992**Date of activities performed:** 9/24/2018-9/30/2018**Date assessment performed:** 9/22/2018**Description of device:** Flexible endoscope**Description of activities performed:** Repair device; lens, irrigation channel, and shaft exterior replaced. Each change was individually and cumulatively assessed.**Determination of whether the activity is remanufacturing:***Lens Assessment*

- Original lens is cracked and needs replacement; OEM lens and epoxy not available for purchase;
- Equivalent lens with same performance specifications and dimensions used (see biocompatibility assessment (BCA) #EI-001 and Component Comparative Analysis Report (CCAR) #EI-002);
- Epoxy used to secure lens is equivalent to OEM epoxy (see BCA #EI-003 and CCAR #EI-004); and
- Leak, optics, and field of view were verified to be within OEM specifications.

*Irrigation Channel Assessment*

- Irrigation channel is worn and leaking fluid into the device;
- OEM part available for purchase and used (part #XX44); and
- Irrigation channel installed and checked for leaks and functionality.

*Shaft Exterior Assessment*

- Shaft exterior damaged during repair activities and needs replacement;
- OEM part not available for purchase; and
- Equivalent shaft exterior with same performance specifications and dimensions used (see BCA #EI-005 and CCAR #EI-006).

*Cumulative Change Assessment*

- Full device specification list inspected and passed (see Customer Evaluation Report #88239 and OEM specification sheet);
- No change in component exposure to reprocessing when following OEM reprocessing instructions;
- A risk-based assessment was performed in each CCAR report; modified risks were identified with using non-OEM parts but were demonstrated as not significantly changing the device's performance or safety specifications, or intended use; and

***Contains Nonbinding Recommendations***

***Draft – Not for Implementation***

- No other change in the risks, or change in the performance or safety specifications, have been identified for the cumulative changes made.

This activity is not remanufacturing.

**Reference to related documents supporting the decision-making process:**

1. BCA #EI-001
2. CCAR #EI-002
3. BCA #EI-003
4. CCAR #EI-004
5. BCA #EI-005
6. CCAR #EI-006
7. Customer Evaluation Report #88239
8. Endoscope Infinity Specification Sheet

**Technician performing service:** xxx

**Reviewed by:** xxx

**Signature(s):** xxx